

Beamline Manual

Version 3.3

S. Stepanyan, E. Pasyuk

August 23, 2023

1 Contacts

Beamline cell phone (757) 303-3996

Personnel listed in Table 1 should be called whenever there is a problem beyond the on-hand expertise. The beam line expert is the main source for help.

Who	Expertise	Cell Phone	Office
MCC	Beam tune, vacuum		x7048 & x7043
Expert on call	Beamline devices and applications	(757) 303-3996	
Eugene Pasyuk	Beamline devices and applications	(757) 515-9811	x6020
Engineering on call	General beamline	(757) 748-5048	

Table 1: The beamline call list.

2 What to Monitor along the Beamline

2.1 Beam conditions

During data taking it is important to monitor electron beam to ensure that there are no changes in the beam parameters which can effect quality of data. For example, changes in the beam halo can dramatically increase background rates making the data worthless, or a change in the electron beam direction can cause unwanted beam losses leading detector damage or beam trip. There are automatic controls and alarms for the beam conditions that will prevent beam damage to equipment and will terminate beam delivery in an event of beam excursions. Nevertheless, the shift taker must continuously monitor the beam parameters and act accordingly.

After the initial beam tune, reference numbers for the beam parameters (position of the beam on position monitors, i.e. BPMs) and the rates of halo counters will be posted in the log book and on the run wiki page, and must be used as a reference.

2.1.1 Beam Current

A consistency in the beam current readings between the Faraday Cup (see Section 3.9) and the "nA" BPMs is one way to check that the electron beam is cleanly transported to the beam dump (Faraday Cup). **Note:** for some of CLAS112 experiments, requested beam current (at nominal luminosity)

exceeds allowed power limit for the Faraday cup, 1 kW (Power in W = Current in nA x Energy in GeV). A beam blocker, see Section 3.9.1, will be inserted in front of the Faraday cup. Even with blocker in-beam, Faraday cup reading is meaningful, it will be a scaled value of the beam current. The coefficient of the proportionality is beam energy dependent [1] and will be calibrated at each given energy setting [2]. With blocker in-beam, the ratio of the currents read by Faraday cup and other current monitors (e.g. BPMs, Section 3.4, and SLM, Section 3.7) must be constant. In all cases one should observe the same current on all three nA BPMs (2C21, 2C24, and 2H01). In addition to the inconsistencies in the beam current readings, the halo counters (see Section 3.3) should show an increased activity if beam scrapes beam pipes. Also higher than normal rates in the detectors would be indicative of scraping immediately downstream of the target.

In the event that there appears to be unacceptable beam losses the following course of action is recommended:

1. Stop beam delivery and data taking, make a log entry flagging any data runs that may have been affected.
2. Call MCC and explain to the operator what has been observed and explain why the tune is unacceptable.
3. Work with the MCC operator to come up with a game plan to fix the problem.
4. Document the solution and start taking data again.

2.1.2 Beam Halo

The presence of a beam halo is usually observed by an increase count rate in the beam halo counters (see Section 3.3). Typically the upstream and midstream counters are very quiet. Any count rate above ~tens of Hz (after initial gain adjustments for a well tuned beam) is indicative of a problem. Note that an increase count rate in the upstream beam halo counters can also indicate an obstacle on the beam path or just a bad beam tune. To further investigate the source of a large count rate in the upstream beam halo counters a harp scan (see the "electron beam profile scan" procedure in Section 3.5) should be performed.

If the beam halo is unacceptable, take the following steps:

1. Stop beam delivery and data taking, make a log entry.
2. Call MCC and explain to the operator what has been observed and why this tune is unacceptable.
3. Work with MCC to solve the problem.
4. Document the solution and start taking data again.

2.1.3 Beam Position

The beam position before CLAS12 is available from three "nA" cavities, 2C21, 2C24, and 2H01, and one stripline BPM, 2H00 (when beam current is above 25 nA). The position monitors close to the target, 2H00 and 2H01, are the most important ones to keep the beam position steady on the target. A feedback system (the orbit locks) uses (x,y) positions on these BPMs, and the Horizontal and Vertical correctors on 2C22/2C23/2H00 girders to keep the beam position stable. Drifts more than ± 0.2 mm should be brought to the attention of MCC. The nA BPMs measure beam position as well as the current. Before data taking, shift worker must confirm that the feedback system is active (verify with MCC). Note, that at beam currents below 25 nA only "nA" cavities are reliable.

2.2 Beamline Vacuum

The beamline vacuum on Hall-B line (from upstream tunnel to the Faraday cup) can be monitored from the vacuum gauge readings available on the Main scaler GUI (see below). The vacuum upstream of the tagger magnet (upstream tunnel) usually is better than 10^{-6} torr, vacuum downstream of the tagger magnet and in the downstream beam line it is usually better than 10^{-5} torr. Vacuum is tightly monitored and interlocked to the beam delivery.

2.2.1 Catastrophic Loss of Vacuum

There are two thin windows that are components of the CLAS beamline, the target scattering chamber exit window and the entrance window to the downstream beamline (after the target). If the window on the scattering chamber fails, fast valves interlocked to the pressure gauges will close automatically. These valves will limit the loss of the vacuum to a small region of the Hall B beamline. Valves are interlocked to the beam Fast Shutdown System (FSD) and beam will be terminated in an event of a vacuum loss.

If any of the valves close due to poor vacuum:

1. notify MCC immediately, turn off the beam (if it is not already OFF)
2. call the engineering on call

If the window on the downstream beamline will fail, rates on the downstream halo counters will increase and will initiate FSD.

2.3 Magnet Power Supplies

The Hall-B beamline and CLAS12 have magnets for beam transport, and momentum analysis. In Table 2 the list of the magnets, their power supplies and point of control (POC) are shown. The items listed with B as the point of control will be controlled by staff in the counting house. The vertical and horizontal correctors, and the tagger magnet are controlled by MCC, but the shift taker should monitor their settings. The tagger magnet power supply is interlocked to the machine Fast Shutdown System (FSD). Interlocks are activated if run requires dumping the beam in the tagger beam dump, e.g during the initial beam tune or during the beam polarization measurements using the Möller polarimeter. If the tagger dipole power supply will trip, the beam will be automatically shut off. This interlock must be masked out when electron beam is delivered to the CLAS target. This is controlled in the BTA GUI by proper state of beam delivery, "photon" or "electron").

The CLAS12 solenoid is also interlocked to the machine Fast Shutdown System and must be activated (unmasked) when electron beam is delivered to the CLAS12 target, to the electron dump (Faraday cup).

If the tagger or any of CLAS12 superconducting magnets do trip off or is set incorrectly take the following action:

1. Call MCC immediately, tell them to shut off the beam (if it is not off already). Call Hall-B engineering on-call if torus or solenoid tripped.
2. Make a log entry.
3. Restore magnets
 - have MCC restore the tagger dipole
 - or

magnet	power supply	POC	Function	Status
MB2C21V		MCC	vertical kick	Active
MB2C21H		MCC	horizontal kick	Active
møller A	Dyna-B	B	Møller polarimeter	Not used
møller B	Dyna-C	B	Møller polarimeter	Not used
MB2C22H		MCC	horizontal kick	Active
MB2C23V		MCC	vertical kick	Active
raster_h1	Danfysik	B	first horizontal target raster	4 Not used
raster_v1	Danfysik	B	first vertical target raster	Not used
tagger	Danfysik	MCC	bend beam to tagger dump	Active
MBD2H00H		MCC	horizontal kick	Active
MBD2H00V		MCC	vertical kick	Active
Solenoid	Danfysik	B	Moller shield and momentum analysis in CD	Not used
Torus	Danfysik	B	CLAS12 FD spectrometer magnet	Active

Table 2: List of the magnets along the Hall B beamline and their functionality.

- call engineering on-call to restore CLAS12 superconducting magnets (torus or solenoid)
4. Restore beam, (if beam was going to Faraday cup, verify that the beam is incident on the downstream viewer on the same location as before the trip.)

2.4 Magnet and Power Supply Beacons

Near every magnet in Hall B there is a **red flashing beacon** that indicates status of the magnets. If beacon is flashing then the magnet **is powered** or **can be powered** at any time. If you need to work near or on a magnet and the red light is flashing you must turn off the supply. The dangers of working near a magnet are limited to those associated with stray magnetic fields. All the high current bus work is enclosed in protective shields so there is no shock hazard. Of course the supply needs to be shut off, and locked and tagged before any of the protective shield is removed.

After any work that required the power supply to be locked and tagged, a through sweep of the magnet area for magnetic debris is required before the lock and tag can be removed.

3 Hall B Beamline devices and Epics Control Screens

All the GUIs necessary to control and monitor beam line devices and the CLAS12 detector can be launched from CS Studio EPICS interface.

3.1 EPICS GUI Launcher

The "*clascss*" (see Fig.1) is a *css studio* screen that serves as an icon manager of *css* GUIs. From this screen, the bulk of the epics applications can be started. It is recommended that when finished with a particular *css* screen, that application be terminated (**not iconized**). If the screen is needed in the future, just launch it again from the "*clascss*". In this manner, searching through all the degenerate icons is eliminated.

To start the "*clascss*":

- log on one of **clon** PCs or **clonsl(1,2,3)** computers as **clasrun**
- type: **classcss**

Most of the beamline applications are under the "*Beamline*" pull down menu. The beamline devices that have motor controls can be found under "*Motors*" pull down menu, scalars under "*Scalars*", the Møller polarimeter controls under *Moeller*.

3.2 The Main Control and Monitoring Screen

The most relevant information for running the beam in a safe and efficient manner is displayed and can be controlled from the Main Control and Monitoring GUI, see Figure 2. This screen can be launch from CS-studio interface using "*Beamline*" button from the menu. The main control screen should be up on a single display on one of PCs in the counting room and must be monitored by shift personnel at all times.

The GUI consists of several rows of information. On the top row rates in the beam halo counters and in the CLAS12 detector elements are displayed. The second row displays information on the beam position (BPM, BOM) and current (BPM, SLM, Faraday cup) monitors, and on the CLAS12 target. The middle row has all moving devices, harps, the collimator, beam viewer, and beam blocker. In row four, vacuum gauges (cold cathode gauges) are displayed. The bottom row shows magnet settings, including correctors, quadrupoles, the tagger dipole and the CLAS12 solenoid and torus. In addition, there are information about the beam, its destination, energy, on helicity related parameters, and the beam viewer in the downstream tunnel.

3.3 Beam Halo Counters (**classc1/classc4**)

The beam halo counters consist of photomultiplier tubes with plastic attached to their photocathodes strapped to the beam pipes along the beamline. There are two halo counters upstream of the Hall-B tagger magnet, *Upstream-L* and *Upstream-R*, two on the top and bottom of the tagger magnet, *Tagger Top* and *Tagger Bottom*, one downstream of the collimator box, *2H00 Down*, four counters upstream of the CLAS target, *Midstream L/R/T/B*, and four counters located in the apex of the forward carriage, *Downstream L/R/T/B*. The beam halo counter scalars are displayed on the Main GUI, top line. HV control of these PMTs can be opened from "*Beamline*" menu.

The rates in halo counters are good indicators of the beam transport quality. After initial beam tune and calibration, rates for good beam conditions will be recorded in the logbook and on the run wiki, and should be closely monitored.

3.3.1 Beam Offset Monitor

In addition to the halo counters, there is a beam offset monitor (BOM) mounted upstream of the target, centered on the entrance window of the target cell. The BOM consists of a glass-silica cylinder (10 mm ID) with 16 plastic fibers attached to the upstream perimeter. A 2 meters long fibers are uniformly distributed along the perimeter of the cylinder, routed outside of the beam vacuum and read out with multi-anode (16) Hamamatsu PMT.

The main function of the BOM is to alert on beam motion at the target. If beam moves from the target center it will generate high rates in BOM PMT due to interaction of the beam tails with quartz cylinder.


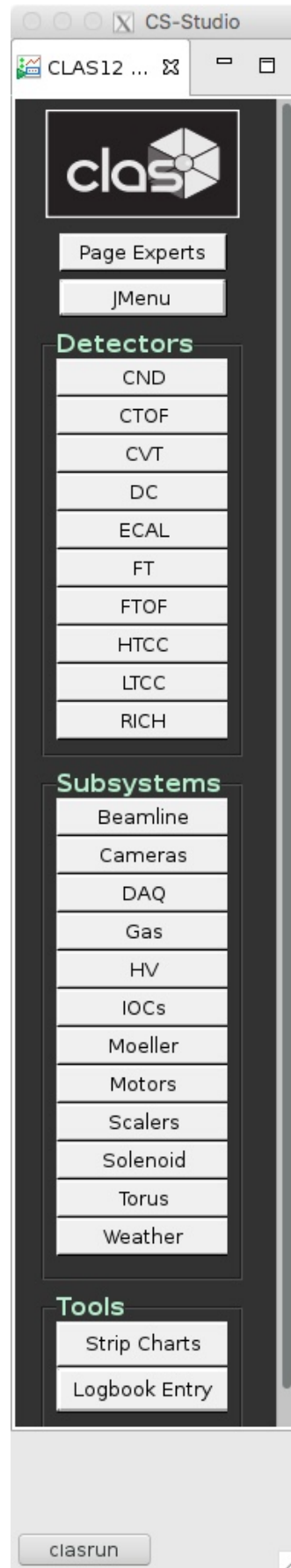


Figure 1: CLAS12 EPICS GUI launcher.

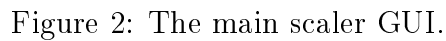

3.4 BPMs

There are 3 nA Beam Position Monitors (BPMs), 2C21, 2C24 and 2H01, and one stripline BPM 2H00 on the B-line. The nA BPMs measure beam transfer (x,y) position relative to their center as well as the beam current. The stripline BPM measures beam position and is reliable only for beam currents above 25 nA. There is a BPM alarm GUI, Fig. 3, that allows to set limits for beam position and current drifts from the requested values. The GUI can be open from the "*Beamline*" pull down menu.

3.5 Hall-B Harps (classc3/classc4)

There are three wire harps on the Hall B beamline, 2C21, 2C24 ("tagger"), and 2H01A. The 2C21 harp has two 25 μm tungsten wires that cross the beam in horizontal (X) and vertical (Y) directions. The "tagger" and 2H01A harps have two sets of wires, thin and thick. There are three thin 25 μm tungsten wires, which cross the beam in X, Y, and 45° axis. This are the main set of wires on these two harps that are used to perform beam profile measurements. Two thick wires (1 mm iron on 2C24 and 0.5 mm carbon on 2H01A) that are mounted farther from the beam, move into beam in X and Y directions and are used only for special measurements of the beam tails. The harp launch GUI, *Harp Scans* shown in Figure 4, enables operator to open controls for the desired harp. Position of the stepper motor that moves the wires in conjunction with the beam halo scalers is used to perform the beam profile measurements. The harp operation is controlled from the *Harp Scans* GUI (controls for 2C21 harp are shown on the figure). During the scan the scaler readout is controlled by the scan application, (see Section 3.3).

In order to perform a harp scan one should push *SCAN* button. After the scan is finished (motor position has been restored at 0.0) use *Analyze* button to see the beam profile and fit the halo counter rate as a function of the wire position. Make sure beam is stable during the scan. If beam has been lost during the scan (e.g. due to beam trips), information (halo counter rates) from wire crossing the beam will be lost. In this case beam scan must be repeated.

3.6 Collimator

In order to prevent direct beam exposure of the SVT or forward MM tracker in an event of beam excursions, a collimator is installed ~ 15 meters upstream of the detectors. This is the same collimator box and a 20 cm long Nickel collimators that have been used in the past for the tagged photon beams in Hall B. One of collimators has been resized to have 20 mm diameter hole. The collimator control GUI will have sizes and motor positions for each collimator in the box, see Figure 5. In order to position required collimator open the *Motors* GUI, select *Hall-B Colli* tab and click on the collimator button you want to move.

3.7 Synchrotron Light Monitor

Synchrotron light monitor (SLM) is a PMT device that measures amount of synchrotron light generated by electron beam in the magnetic field of the last dipole right before the Hall-B upstream tunnel (when beam gets to the hall beamline elevation level). The synchrotron light device, nomenclature name 2C20, is composed of a mirror and a prism that splits the synchrotron light image in two. One is used by accelerator operations as beam viewer, another is directed to a 2" PMT and read out on current integrator. The PMT signal is used by Hall-B as beam current monitor.

The main use of this current monitor is to measure helicity related beam intensity (charge) asymmetry. While PMT signal strength depends on the beam position, for fast ($>\text{Hz}$) asymmetry measurements

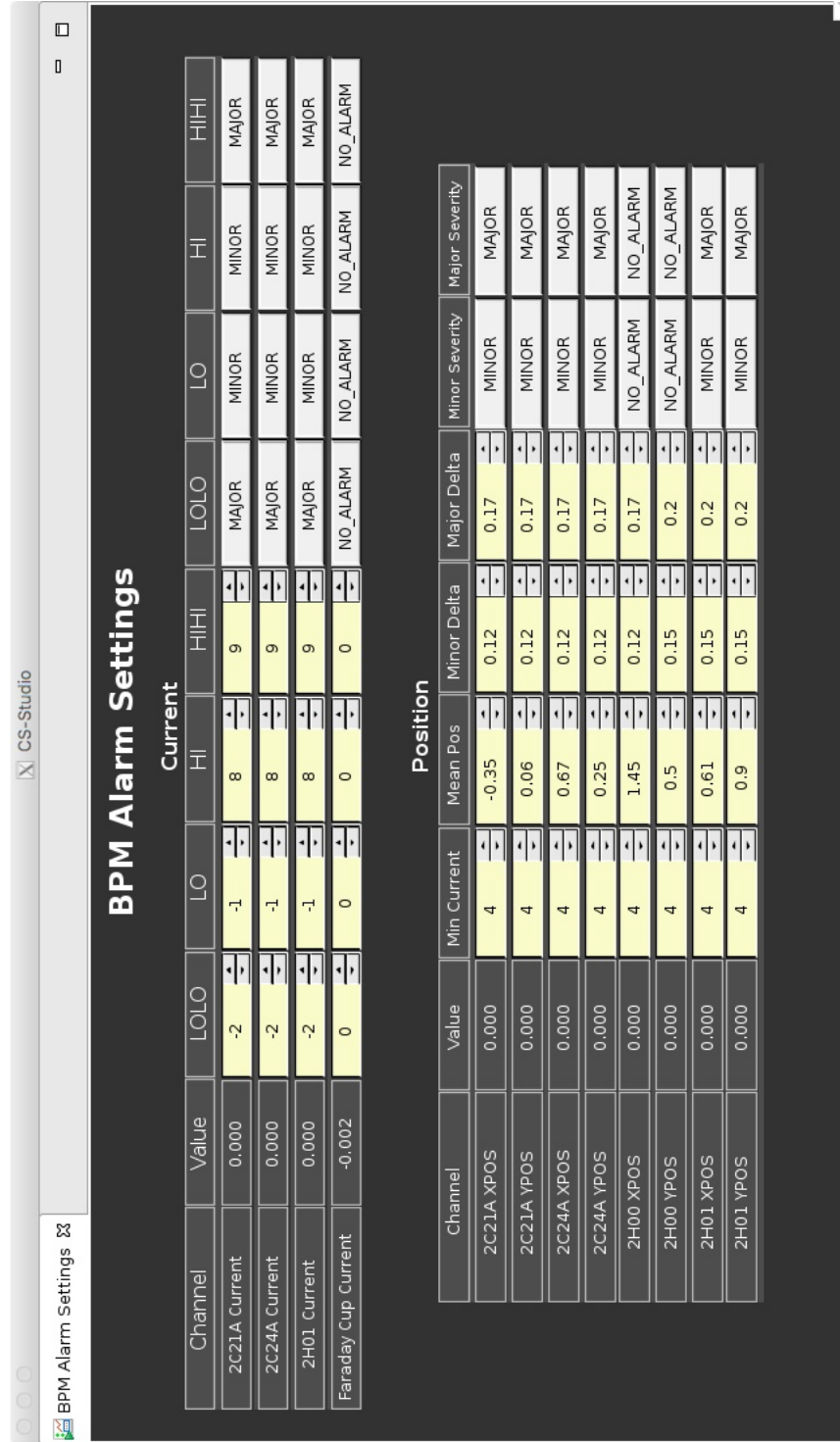


Figure 3: The Hall B BPM setup screen. The current reading, in nAmps. The Faraday Cup current reading is updated at the same rate as the beam halo scalers.

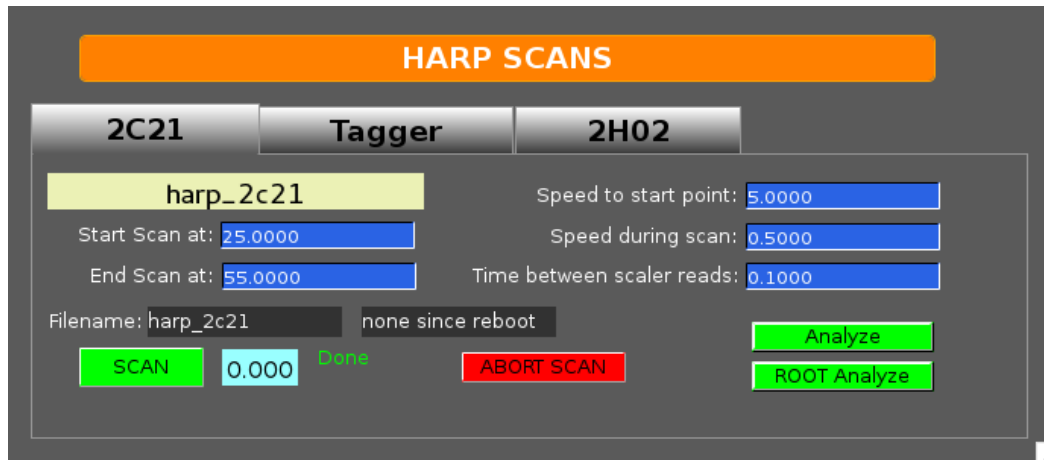


Figure 4: The main Harp medm screen. The green buttons on the top are for opening individual harp control. Cyan buttons below for expert GUIs.

it works well without any corrections for position drifts, since position drifts are usually very slow process. Position dependence can be corrected using the 2C21 BPM position readout and in principal this device can be used for beam charge accounting if needed.

3.8 Beam Viewer

There is a beam viewer installed in the downstream tunnel, right before the shielding wall that separates Faraday cup dump. The viewer is a wheel with a Chromox screen mounted on the axis. The motor that control the position of the screens can be controlled from the same *Motors* GUI as collimator, see Figure 5. In order to position the screen just push the button carrying the name of the screen.

3.9 Faraday Cup (classc4)

The instantaneous beam current reading from the Faraday cup is available on the **Main Scaler** screen as shown in Figure 2. The update rate is the same as for the beam halo scalers and is controlled from the scaler GUI, see Section 3.3. An important consideration is that the Faraday Cup current integrator rate is 906.2 counts/sec for 1 nA. This means that if the count time on the scaler is much less than 1 sec you will observe large statistical fluctuations.

3.9.1 Beam blocker

The Hall-B Faraday cup is not cooled and cannot operate at high currents (deposited hit should not exceed 1 kW for more than 1/2 hour running). If run requires use of beam currents above the limit (for 11 GeV it is 90 nA), a 17 kW beam blocker, cooled-copper absorber, must be position in front of the Faraday cup. The beam blocker can be controlled from *Motors* GUI shown in Figure 5. From *Motors* GUI select *Blocker* then push *Go beam* button in order to put the blocker on the beam and *Go Home* to retract if from the beamline.

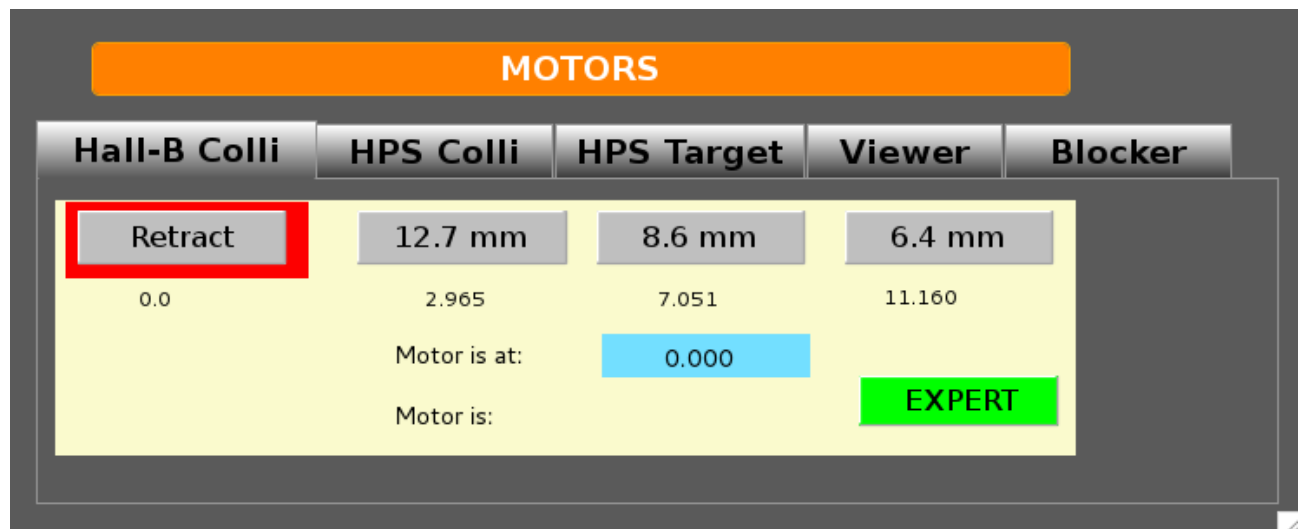


Figure 5: The beam GUI to operate collimator and beam blocker.

3.10 Target

The Hall-B cryo target is used for the CLAS12. The target cell diameter is designed to be 20 mm and the thin portion of the cell windows (30 μm thick aluminum) 10 mm in diameter. Figure 6 shows design rendering of the CLAS cryo target inside the foam scattering chamber. Control of the target, filling and emptying, will be done using *Target PC*. Operator should follow special instructions for the target operations. For monitoring, the target temperature, pressure, and the buffer dewar liquid level are displayed on the Main Control and Monitoring GUI, Figure 2.

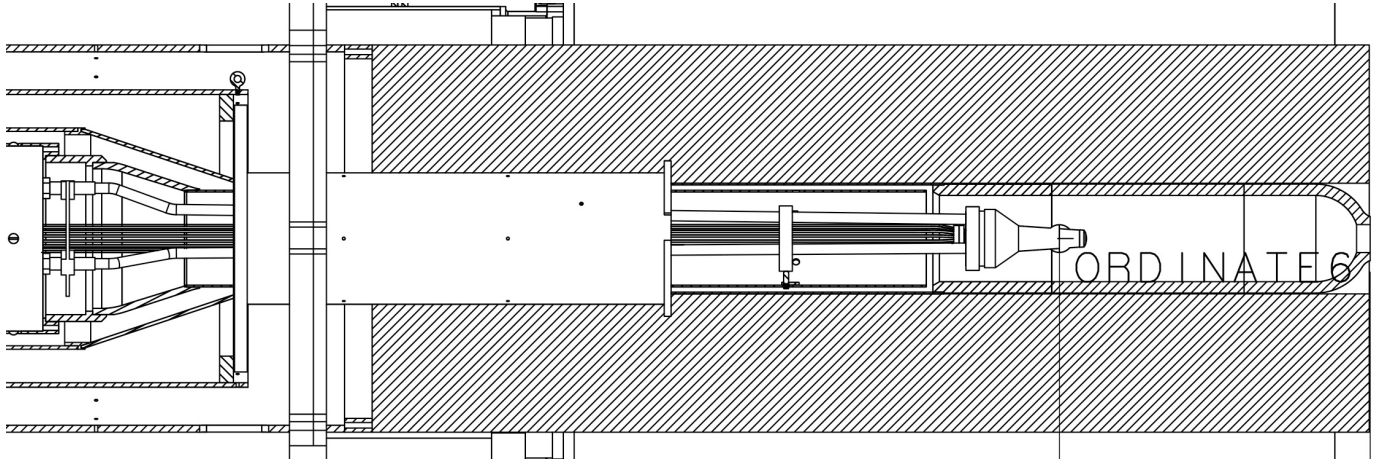


Figure 6: Hall B cryo target and scattering chamber.

Appendix

4 Establishing beam for physics

This section describes a generic beam tuning procedure. Each experiment may have special requirements which are described in the experiment specific procedure.

Establishing production quality electron beam for experiments in Hall B is a two step process.

4.1 Beam to tagger (tagger-yoke) dump

The initial tune is done at low currents, < 10 nA. If the beam energy is below 6.12 GeV it is dumped on the tagger dump located under the hall's floor. If the beam energy is above 6.12 GeV it is dumped on the "tagger-yoke dump", a dump in the tagger dipole magnet yoke, before hall proper. Beam is deflected down to this intermediate dump by the tagger dipole magnet [3]. The tagger dipole power supply set current relates to the beam energy as [4]:

$$I(A) = 43.491 \times E(\text{GeV}) - 0.076 \quad (1)$$

In this first step:

- the CLAS12 detectors (especially tracking detectors) must be OFF
- all halo counters must be ON
- masked the halo counters (upstream, midstream, BOM and downstream) in beam Fast Shut Down (FSD) system
- the "blank" collimator with no hole is on the beam
- the CLAS12 solenoid and torus magnets are energized (or can be energized while beam tune is in progress)

The beam with required profile and trajectory is established by MCC OPS using correctors and quadrupoles on 2C line, while monitored by Hall-B shift crew using the wire harps and the nanoamp (nA) BPMs [5] at 2C21 and 2C24 girders in the upstream tunnel (see Figure 7). There is a Yag viewer, ITV2C24, upstream of the tagger dipole, controlled by MCC, that may be used by MCC to verify position and profile of the beam. If the beam is tuned to the tagger dump the beam spot can be seen on the tagger dump viewer (see Figure 8). Hall-B shift personnel should work with MCC operator to establish required quality beam. All harp scans must be properly analyzed and logged into logbook. The beam tune is good when required parameters for the profile (x/y widths) and trajectory (x/y positions on various monitors) are achieved. These parameters will be written on the run wiki and/or on the white board in the counting house. Typically widths at 2C21 harp $\sigma_{x,y} \leq 200 \mu\text{m}$, while at 2C24 (tagger) harp $\sigma_{x,y} \leq 500 - 600 \mu\text{m}$. In general it is a good practice to check newly measured beam parameters against last good/acceptable tune.

Description	Name	L (meters)
SLM	ISR2C20	-44.3
Green shield wall		
Stripline BPM	IPM2C21	-41.4
Quadrupole	MQR2C21	-41.1
Horizontal corrector	MBC2C21H	-40.6
Vertical corrector	MBC2C21V	-40.1
Quadrupole	MQA2C21A	-39.6
Horizontal corrector	MBC2C21AH	-39.1
Beam viewer	ITV2C21	-38.9
Wire Harp	IHA2C21	-38.8
nA-BPM	IPM2C21A	-37.6
Stripline BPM	IPM2C22	-26.9
Quadrupole	MQK2C22	-26.5
Horizontal corrector	MBC2C22H	-26.2
Quadrupole	MQK2C23	-25.8
Vertical corrector	MBC2C23V	-25.5
Quadrupole	MQK2C24	-24.9
nA-BPM	IPM2C24A	-24.5
Wire Harp	IHA2C24	-22.0
Beam viewer	ITV2C24	-21.8
Hall-B tagger dipole	TAGGERB	-17.6
Hall-B collimator	ETA2H00	-17.0
Stripline BPM	IPM2H00	-12.3
Quadrupole	MQA2H00	-11.9
Quadrupole	MQA2H00A	-11.6
Horizontal corrector	MBD2H00H	-11.3
Vertical corrector	MBD2H00V	-11.1
Quadrupole	MQB2H01	-8.6
nA-BPM	IPM2H01	-8.0
Wire harp	IHA2H01A	-7.5
Center of the hall		0
Stripline BPM	IPM2H02	13.5
SVT collimator	ETA2H02	14.1
Wire harp	IHA2H02A	14.8
Dipole 1	MFC2H02A	15.3
HPS target	ETA2HHPS	17.0
Spectrometer Dipole	MFC2H02B	17.5
Dipole 2	MFC2H02C	19.7
Beam viewer	ITV2H04	24.0
Dump, Faraday cup	IFY2H04	27.0

Figure 7: Bernaline elements from the green shield wall to Faraday cup dump.

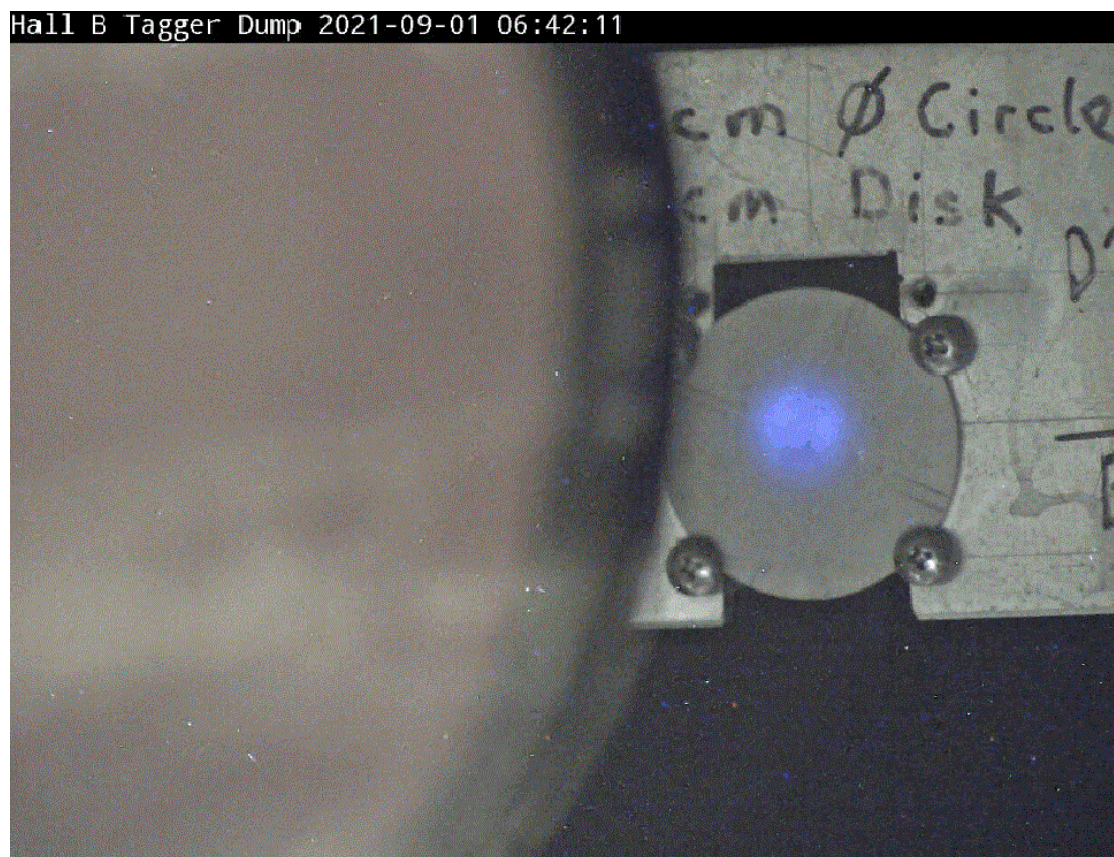


Figure 8: Beam spot on the tagger dump viewer.

4.2 Beam to Faraday cup

The second step of establishing a physics quality beam starts after acceptable beam parameters have been achieved on the tagger-yoke dump. The following are steps for sending beam to CLAS12:

- tell MCC that beam is acceptable at the tagger
- ask to take the beam away, and degauss and turn OFF the tagger magnet
- position 20 mm diameter collimator on the beam
- set the halo counter FSD thresholds in to 1 MHz and the integration time interval to 50 milliseconds
- when ready ask MCC to send ~ 5 nA beam straight to the electron dump, at the end of the Hall B beamline where Faraday cup is located
- verify that beam goes to the dump
 - a. make sure beam is clearly visible on the downstream viewer, use Chromox screen¹. Beam should be within 10 mm (2-tick marks) around the center (see Figure 9)
 - b. make sure that Faraday cup beam current reading and the beam currents on 2C21 and 2C24 BPMs are consistent (the difference should not be more then few %)

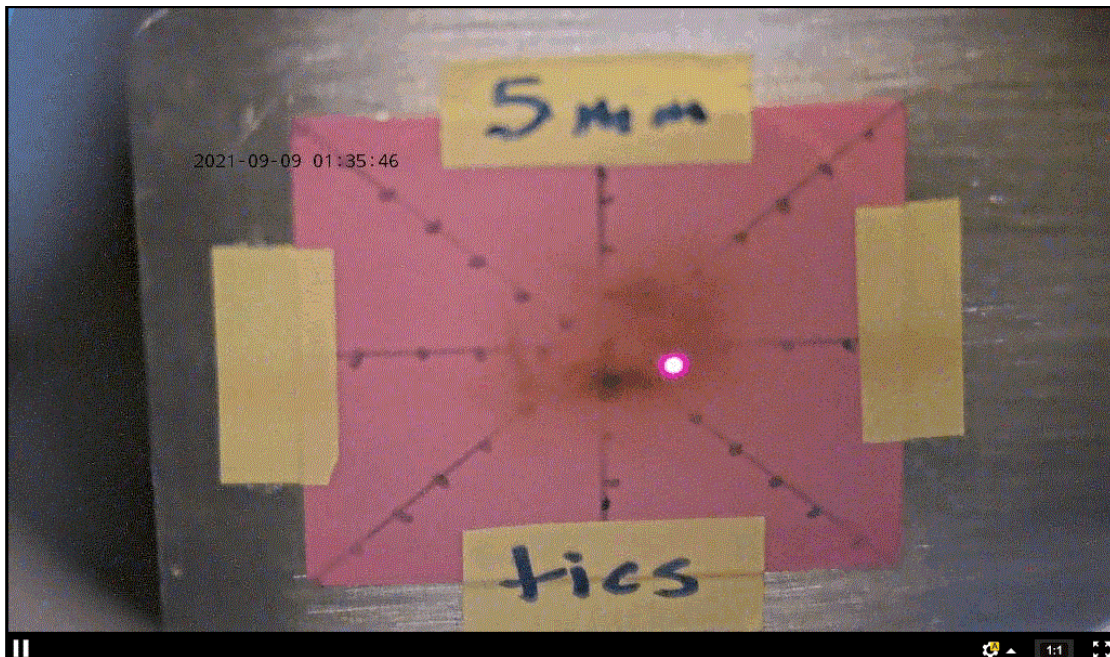


Figure 9: Beam spot on the Faraday cup viewer.

The beam profile and position adjustments on the target will be done using correctors and quadrupoles on 2C22/2C23/2C24 girders in the upstream tunnel and 2H00 girder in the hall (if needed). The last one is the closest to the target (~ 10 m upstream) and will be used to focus beam at the target location to achieve required size (preferably $< 300 \mu\text{m}$). The profile and position of the beam on the CLAS12 target will be checked using a 3-wire harp 2H01A mounted about ~ 7 meters upstream of the CLAS12 target and 2H01 nA BPM (stripline BPM on 2H00 will not work at low currents). The 2H01A harp

¹At low energies beam spot may not be clearly visible at 5 nA

measures the beam profile and its projected position along x-, y-, and 45° axes. After physics quality beam is established based on the profile, beam position on either the cryo target cell or any of CLAS12 detector component (e.g. FTcal) must be adjusted based on the lowest rates on the downstream halo counters and on BOM or symmetry of rates in the detector. Procedure for adjusting the beam on the target is the following:

- ask MCC operator to move beam on 2H01 BPM by 0.1 mm steps up and down, vertical direction
- record rates on halo counters and on BOM for each step. Stop moving in the given direction when rates go more more than $\times 10$ from where you started
- analyze rates as a function of position, find a position on 2H01 corresponding to the mid point of the two extreme ends where rates were highest. Ask MCC operator to position the beam on that position on 2H01 nA BPM
- repeat everything for horizontal alignment, moving the beam to left and right
- set orbit lock using found vertical and horizontal positions on 2H01

For adjusting the beam position to get symmetric rates on a detector, e.g. FTcal, move beam with small steps in x- and/or y- on 2H01 nA BPM until desired symmetry is achieved.

After high quality beam was established and properly aligned, the beam orbit lock system should be engaged. This system incorporates position readings from the BPMs (2H01 and 2H00 if currents are high enough) to regulate currents in the horizontal and vertical corrector dipoles to minimize beam motion at the target. The final step in establishing the production running conditions is setting limits on the halo counter and BOM rates for FSD system. If beam moves unexpectedly and gets close to an obstacle, e.g. collimator walls or to the thick parts of the target cell, count rates on the beam halo monitors and on the BOM will increase. The appropriate rate limits will depend on actual run conditions and the target, and will be noted on the white board in the counting room and/or on the run wiki. General prescription for setting a limit for FSD input rate of $N(\text{Hz})$ for a trip time interval of δt is:

$$N_{thr} = N + n_{\sigma} \times \sqrt{\frac{N}{\delta t}}$$

where n_{σ} is how far from mean value we want the trip to acquire. For example, for $\delta t = 5$ ms - the recommended value, $n_{\sigma} = 5$. will mean one false tripe every 5 hour. In order to avoid frequent false trips and allow some overhead in average rate and possible smaller time interval, a recommended value is $n_{\sigma} = 6$. Note, when beam is cleanly transported to Faraday cup, "Upstream" and "Midstream" halo counters should count $< 10 - 15$ Hz, for those FSD threshold should be set to 1000. High rates (~ 100 Hz) on these counters will indicate bed beam transport or a bleed-through from other halls and must be corrected.

5 EXPERT procedure for beam polarization measurement using the Möller polarimeter

This section is for experts. For routine operation by shift operator follow the procedure described in "Performing Möller Runs" document.

5.1 Before start

- Reasonable quality beam must be already present for Hall-B
- Beam has to be terminated on the tagger-yoke dump
- CLAS12 is OFF, especially sensitive detectors: HV on **Drift Chambers** and **SVT/MVT**

5.2 Setup

1. Notify MCC that you are about to do Möller run and **request to take the beam to the tagger yoke beam dump** (they will need to take the beam away and energize the tagger magnet), MCC will ask to change BTA setting to "photon"
2. Ask MCC to turn orbit locks off, and mask BOM and Halo Counters in FSD
3. Turning ON the polarimeter is done from EPICS GUIs (for now multiple control GUIs in expert mode are used).

NOTE: BEAM SHOULD BE OFF DURING the SETUP of MöLLER POLARIMETER
Möller setup GUIs can be launched from the **Moeller** tab on the *"clascss"* GUI.

- (a) the "Moeller Asym - All" GUI, see Fig.10 contains all the helicity gated scalers, charge asymmetry, and the beam polarization². It has several controls for the measurement and monitoring:

Buttons "Start", "Reset" and "Stop" will start and stop acquisition of data or reset acquisition (clear scaler buffers).

The acquisition time controls update frequency of scalers, measured polarization value, and the charge asymmetry value. It is recommended to have Möller polarimeter acquisition time > 60 seconds when taking the measurements. For practical reasons, at the start when charge asymmetry will need tuning (see below) time for the scaler update can be set to ~ 10 seconds. The time can be set either by typing a value in the box or moving the slider.

The "Attenuator Controls" are to control beam charge asymmetry. By changing voltage on intensity attenuators (IA) one can equalizes intensity across helicity states. It is recommended to use "Global Offset" that will change voltage on all four IAs at the same amount. Desire charge asymmetry is $< 0.1\%$ based on SLM.

- (b) The control for Möller quadrupole power supplies are provided in "Moeller Quadrupoles" GUI, see Fig.11. Power supplies will be turned ON and in remote before hall closing. From GUI one should first turn them on by pushing the "PS ON" buttons, then set the desire value for currents in "Current Setpoint" window. For 10.7 GeV the suggested value for the quadrupoles is 3050 A.
- (c) The target is polarized to its saturation by a longitudinal (along the beam) magnetic field generated using pair of Helmholtz coils. It is expected that the target will be saturated at ~ 1.8 A current in the coils. The recommended current for Möller measurement is 3 A. A GUI for power supply of Helmholtz coils, "Moeller Helmholtz PS" see Fig.12, has two controls, button "STATE" defines state of the power supply. Typically it will be in "STANDBY" state when is not used. To energize coils first from the menu in "STATE"

²The polarization and charge asymmetry have GUIs their own, but for now this main GUI will be used.

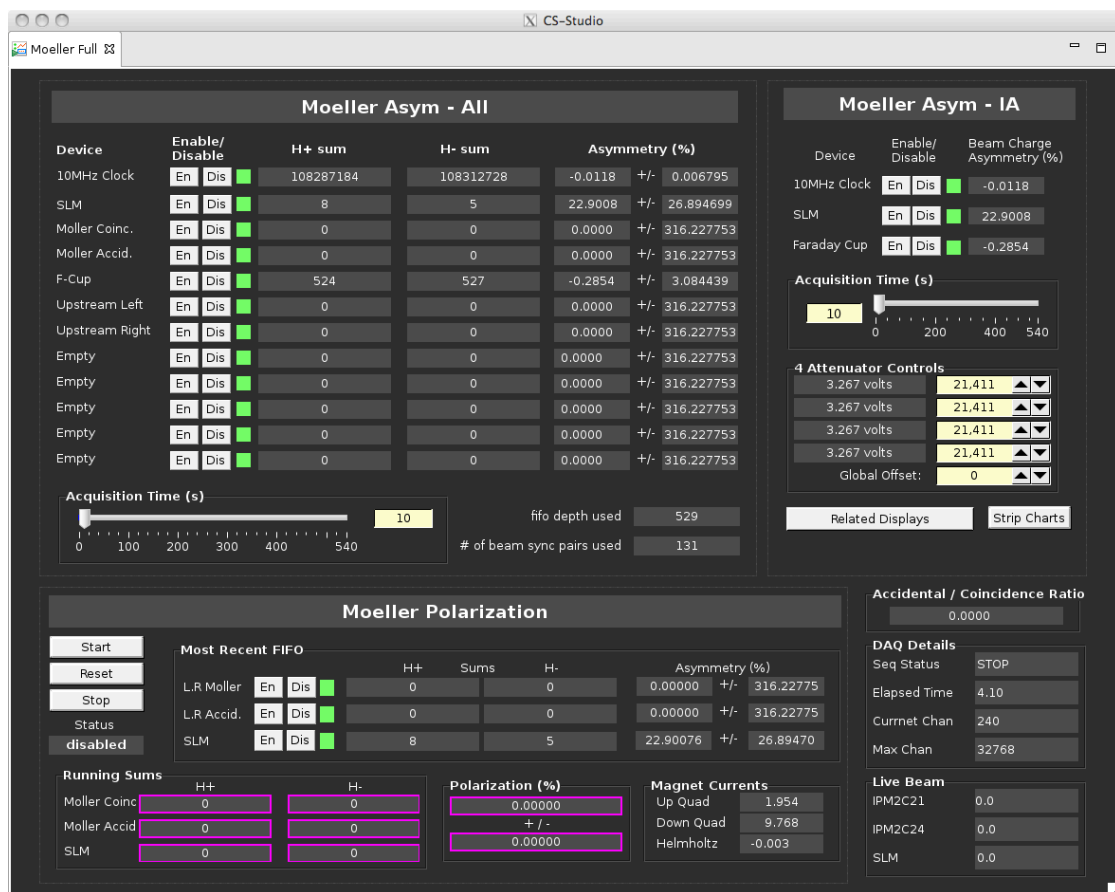


Figure 10: The main Möller Epics expert GUI.

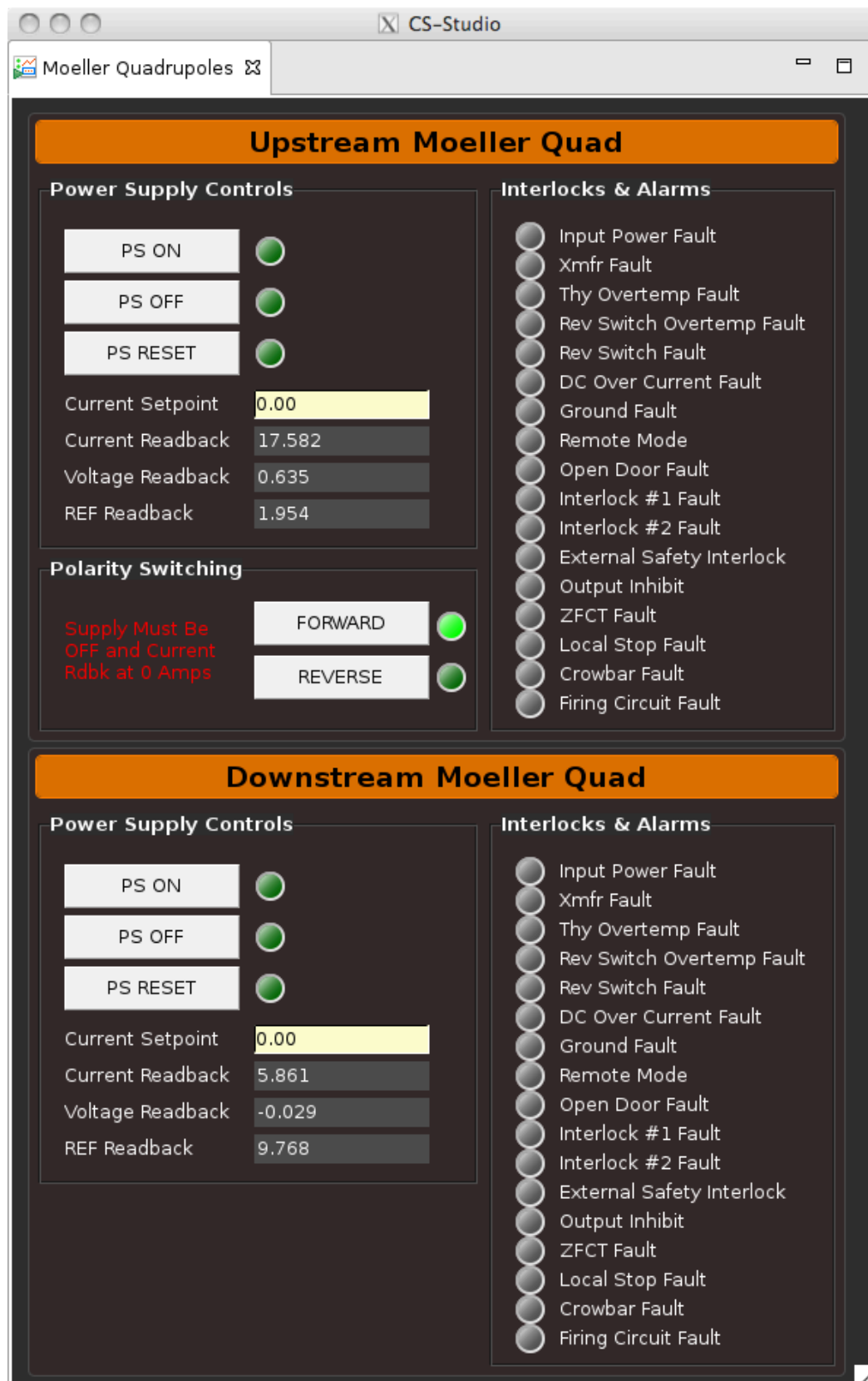


Figure 11: Control GUI for the Möller quadrupole power supplies.

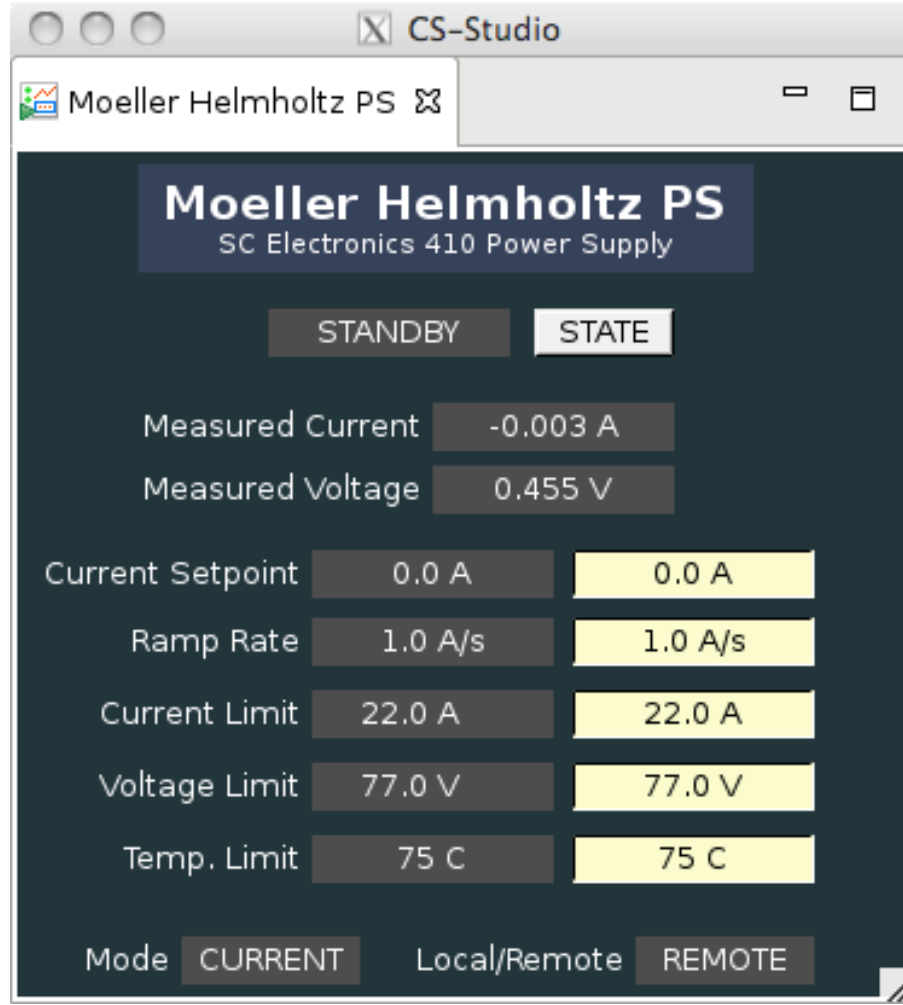


Figure 12: Control GUI for the Möller target Helmholtz coils power supply.

chose PS ON, then in "Current Setpoint", a white window, write the value, either 3 or -3 . Beam polarization measurements with both orientations the Helmholtz field is recommended to check systematics.

- (d) a target control GUI, Fig.13, allows to position desired target on the beam. Left target is the recommended target for the measurements.
- 4. Once the tagger magnet is energized, Møller setup is UP and ready, request the beam current as specify for the given energy Møller measurements to be delivered³, as measured by 2C21 nA BPM and/or SLM. Do not use 2C24A since that BPM is located downstream of the Möller setup. For ~ 11 GeV beam, if beam conditions are normal, as expected, the beam current should be 4 nA.

5.3 Data Taking

³The optimal beam current is a function of beam energy. More specific information may be available on the white board in the counting house or in the run period specific documentation on the run wiki. Regardless of what currents are specified on the white board or in this document, the ratio $Left \otimes Right$ accidentals to the true coincidence rate should be kept below 5%. It may be necessary to adjust the HV on the Left and Right PMT's to achieve a low accidental rate, while maintaining a reasonable true rate.

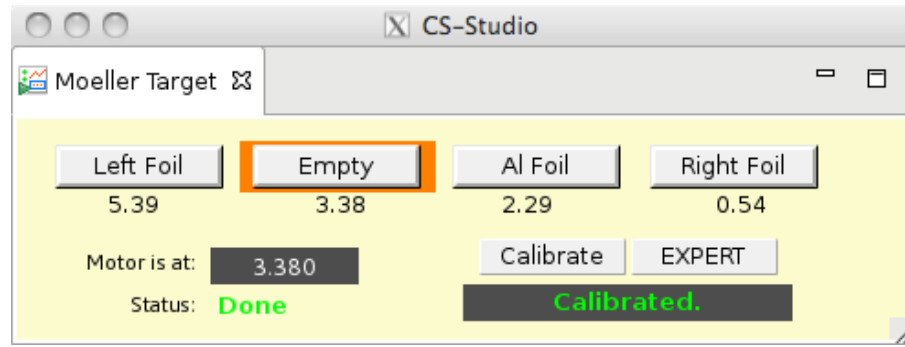


Figure 13: Control GUI for the Möller target.

- To start a new run when DAQ is still running hit the "Reset" button on "Moeller Asym - All" GUI. If run was stopped hit the "Start" then "Reset".
- Run is complete when the error on the beam polarization on the GUI is below $\leq 1.5\%$ absolute. Typically it takes about 45 min to 60 min to get the required accuracy (beam condition dependent).
- Make measurements with both positions of the half-wave plate, "IN" and "OUT" (start the first measurement with whichever position it is, then do the second measurement with the other setting).
- If needed perform measurements with both polarity of the Helmholtz coils.
- Log every measurement by sending "Moeller Asym - All" GUI to logbook together with main scaler GUI to document beam currents, beam position, and halo counter rates.

5.4 Backing off Möller setup out

When done with the measurements:

- Do not forget to make a log entry including all details and the GUI!
- Request MCC to take the beam away and **de-gauss the tagger magnet if the next step is to send the beam to Faraday cup (usually it is).**
- Turn off quadrupoles by setting 0 in "Current Setpoint"s and then when current readback is at ~ 0 A push "PS OFF" button
- Turn off Helmholtz coils by setting 0 in the current "Current Setpoint" and change "STSTE" to "STANDBY" when "Measured Current" is ~ 0
- Retract the target by pushing "Empty" button on the target GUI
- Once tagger magnet is degaussed, restore beam to the Faraday Cup.
- Turn ON CLAS12

References

- [1] N. Dashyan and S. Stepanyan, CLAS12 Note 2016-004
- [2] R. Paremuzyan and S. Stepanyan, CLAS12 Note 2018-003 (2018).
- [3] D. Sober et al., Nucl. Inst. and Meth. A 440, 263 (2000).
- [4] https://clasweb.jlab.org/wiki/images/9/97/Beam_on_tagger.pdf
- [5] M. Piller, et al., JLAB-ACC-99-30 (1998).

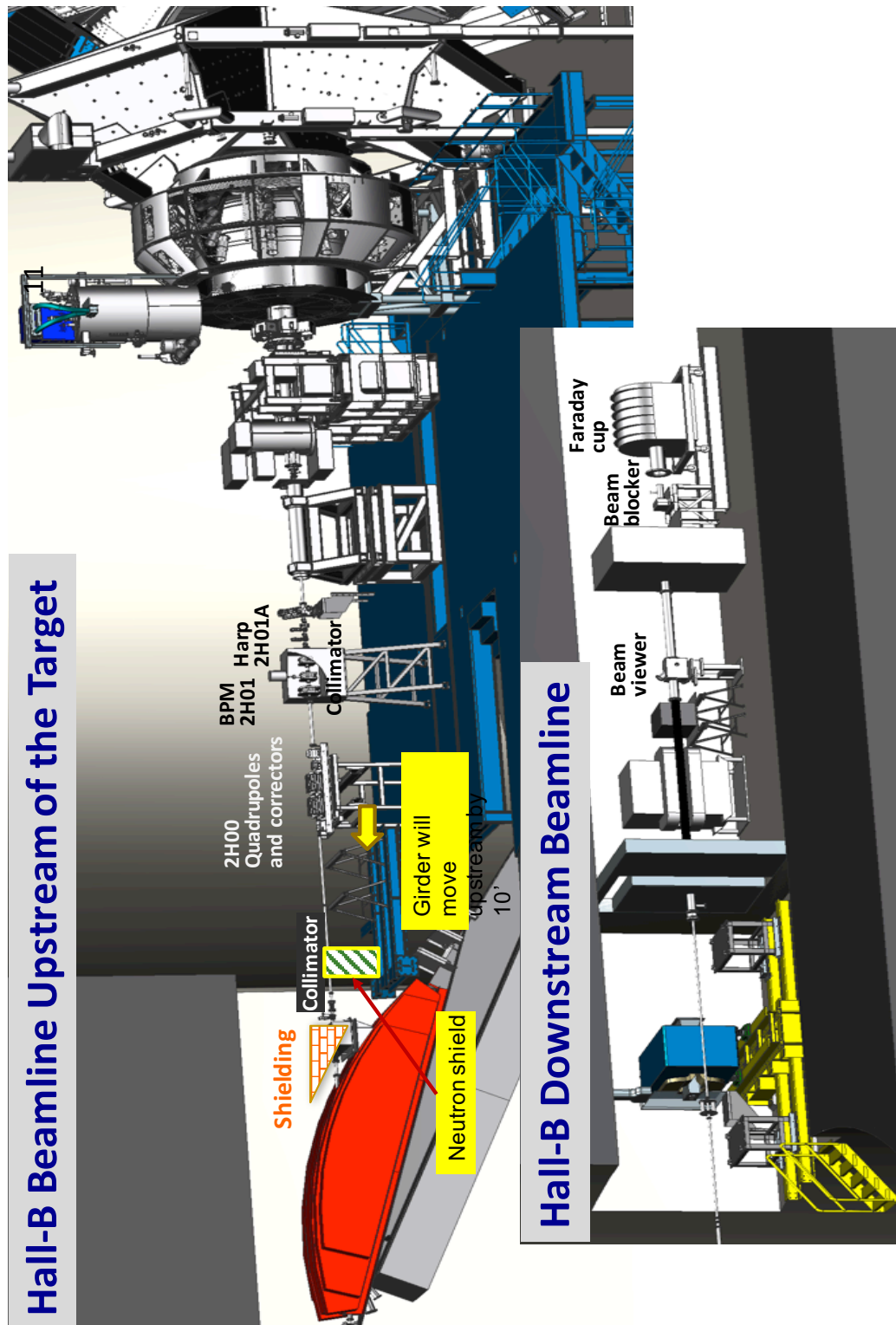


Figure 14: Hall-B beamline.