

CLAS12 Slow Controls Operations Manual - v1.6

(Dated: January 29, 2017)

I. Overview

The operator interface for the Hall B controls systems is based on Control System Studio (also called CS-Studio or CSS) and allows access to all the necessary EPICS tools from a single application. This system is accessible locally by user `clasrun` from all `clonpc##` desktop computers in the Hall B Counting Room for shift workers. For remote access, see Section X.

To start the control system with only the main menu as shown in the left panel of Figure 1, in a terminal run:

```
clascss
```

To start the control system with the full alarm handler as shown in the right panel of Figure 1, in a terminal run:

```
clascss-alarm
```

Before opening CS-Studio on a clonpc desktop in the counting house, check whether an instance is already running on that machine. If one is, use that instance instead of opening another.

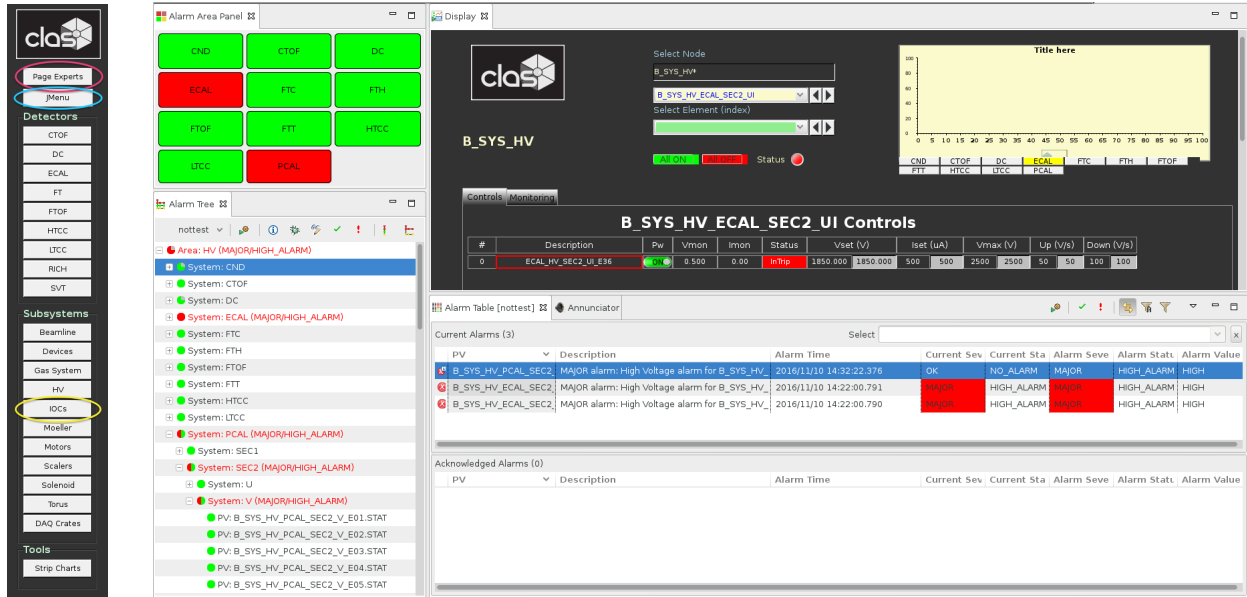


FIG. 1: On the left is the main CLAS12 controls menu showing locations of buttons used for the accelerator's JMenu screens (blue), paging Hall-B experts (pink), and restarting IOCs (yellow). On the right is the alarm controls screen.


The main menu is organized in a hierarchy of subsystems and components. The top portion of the menu is for specific detectors, while the bottom portion is for more general subsystems. The alarm screen is organized into different regions described in the following sections.

II. Alarms

The user frontend of the alarm handling system runs in CS-Studio and includes visual and audible alarms. The default alarm view when running `clascss-alarm` is shown in the right panel of Figure 1 and contains four sections:

1. Top Left: the *Area Panel*, which gives an overview of the global alarm system status. The color of the areas reflects the most severe alarm in that area.
2. Bottom Left: the *Alarm Tree*, a hierarchical view of all alarm settings.
3. Top Right: the *Alarm Table* (see also Figure 2), a list of active alarms that need to be addressed.
4. Bottom Right: the *Acknowledged Alarms* (see also Figure 2), a separate list of active alarms that have already been acknowledged.

When an alarm triggers, it will enter the *Alarm Table* and its color will change according to its severity. The annunciator (running on `clonpc17`) will also audibly announce any new alarms or a count of currently active alarms.

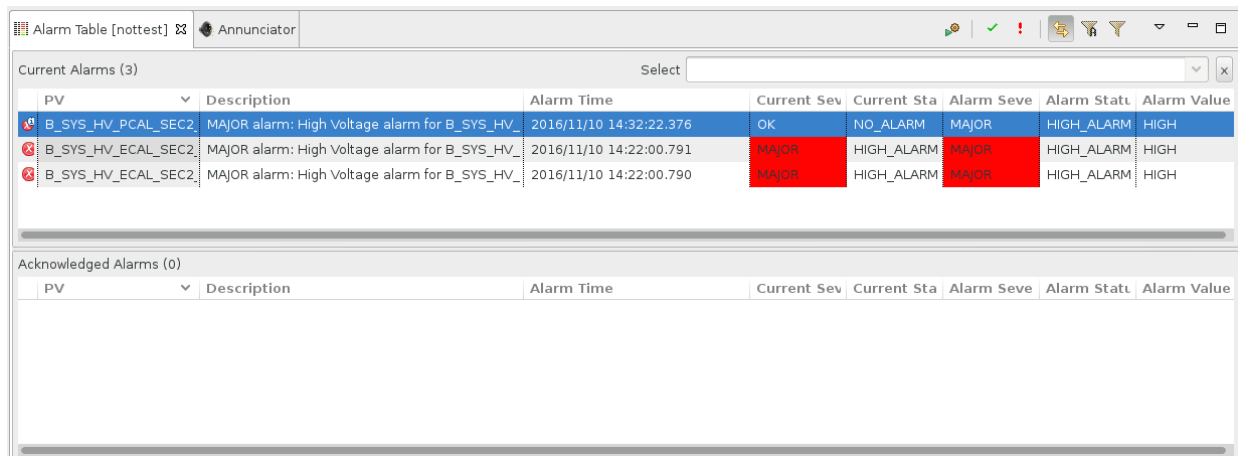

PV	Description	Alarm Time	Current Sev	Current Sta	Alarm Seve	Alarm Statu	Alarm Value
B_SYS_HV_PCAL_SEC2	MAJOR alarm: High Voltage alarm for B_SYS_HV_	2016/11/10 14:32:22.376	OK	NO_ALARM	MAJOR	HIGH_ALARM	HIGH
B_SYS_HV_ECAL_SEC2	MAJOR alarm: High Voltage alarm for B_SYS_HV_	2016/11/10 14:22:00.791	MAJOR	HIGH_ALARM	MAJOR	HIGH_ALARM	HIGH
B_SYS_HV_ECAL_SEC2	MAJOR alarm: High Voltage alarm for B_SYS_HV_	2016/11/10 14:22:00.790	MAJOR	HIGH_ALARM	MAJOR	HIGH_ALARM	HIGH

PV	Description	Alarm Time	Current Sev	Current Sta	Alarm Seve	Alarm Statu	Alarm Value
----	-------------	------------	-------------	-------------	------------	-------------	-------------

FIG. 2: The *Alarm Tree* portion of the alarm screen, showing an example of three outstanding alarms that need to be addressed. The first is no longer in an alarm state (denoted by the *OK* in the *Current Severity* column), and none of the three have been acknowledged (or else they would have appeared instead in the lower *Acknowledged Alarms* section).

By right-clicking on an alarm in the *Alarm Table*, a dropdown menu of actions is accessible (see Figure 3). This dropdown list contains access to a *Guidance* screen with instructions that should be read and followed on how to deal with the specific alarm. The next step is to acknowledge the alarm using the *Acknowledge* option in the dropdown menu, which will silence the alarm and move it to the *Acknowledged Alarms* section until it is no longer in an alarm state. For many alarms there is also an option in the dropdown menu starting with *Open* that will open a screen necessary to address the specific alarm using the information from the *Guidance* screen.

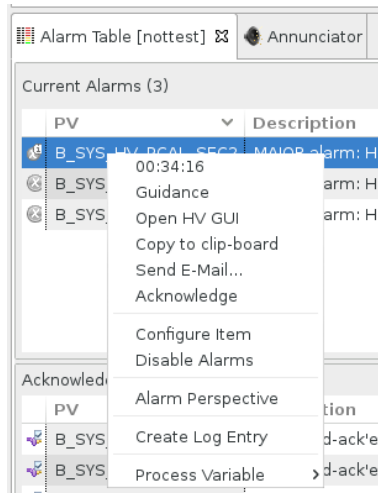


FIG. 3: An example dropdown menu accessible by right-clicking on an alarm in the *Alarm Table*. Important visible actions include a *Guidance* button, an *Open* screen action, and the *Acknowledge* action.

III. IOCs

EPICS input-output controllers (IOCs) are the backend responsible for the actual communication with the hardware devices in the hall. Figure 4 illustrates access to the IOC controls screens from the main CLAS12 menu, as well as the overview IOC heartbeat screen. The heartbeats should be flashing at 1 Hz for all IOCs, or else the IOC may be in need of reboot. By clicking on the IOC in the heartbeat screen (or the IOC health group in the main menu), controls to monitor and reboot the IOCs can be accessed, and an example is shown in Figure 5. Systems are in place to automatically start all necessary IOCs if for any reason they are not running (e.g. recovery from a power outage), however occasionally a manual reboot is required.

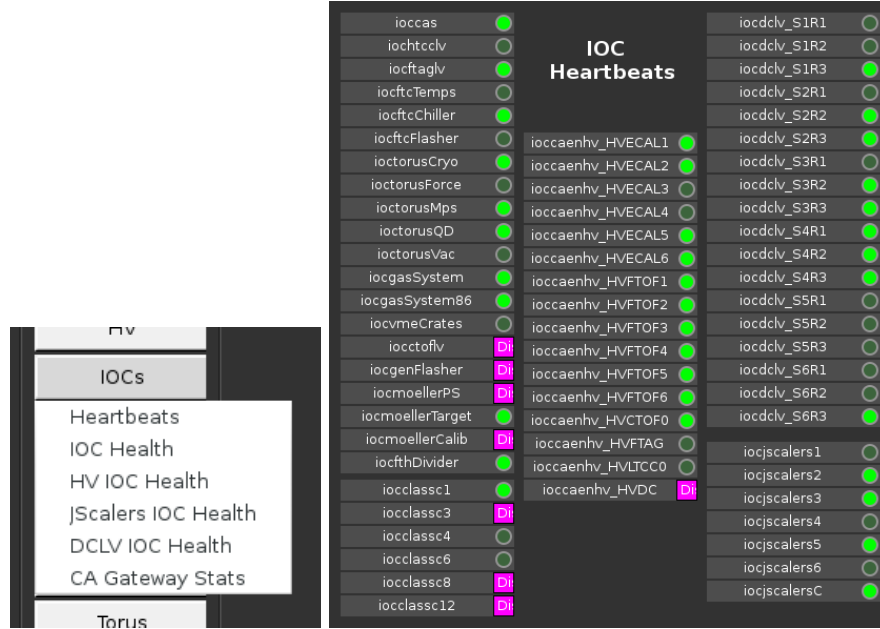


FIG. 4: Dropdown menu (left) from the *IOCs* button in the main CLAS12 controls menu showing links to health screens for subsets of IOC groups, and the IOC heartbeat overview screen (right).

IOC Health											
softIOCs											
IOC Name	Hostname	Up Time	Heartbeat	Expert	Soft Reboot	Last Reboot	Console	Hard Reboot	Status	Message	Autosave
iocccas	clonioc1.jlab.org	15 days, 00:19:00	1297140			10/25/2016 13:57:47			Warning	2 values not saved	Write 'cas_settings.sav1'
iocdclv	clonioc2.jlab.org	43 days, 21:55:55	3794155			09/26/2016 16:20:51			Warning	2 values not saved	Write 'cas_settings.sav1'
iocdclv_S1R1									Warning	2 values not saved	Write 'cas_settings.sav1'
iocdclv_S1R2									Warning	2 values not saved	Write 'cas_settings.sav1'
iocdclv_S1R3									Warning	2 values not saved	Write 'cas_settings.sav1'
iocdclv_S2R1									Warning	2 values not saved	Write 'cas_settings.sav1'
iocdclv_S2R2									Warning	2 values not saved	Write 'cas_settings.sav1'
iocdclv_S2R3									Warning	2 values not saved	Write 'cas_settings.sav1'
iocdclv_S3R1									Warning	2 values not saved	Write 'cas_settings.sav1'
iocdclv_S3R2									Warning	2 values not saved	Write 'cas_settings.sav1'
iocdclv_S3R3									Warning	2 values not saved	Write 'cas_settings.sav1'
iocdclv_S4R1									Warning	2 values not saved	Write 'cas_settings.sav1'
iocdclv_S4R2									Warning	2 values not saved	Write 'cas_settings.sav1'
iocdclv_S4R3									Warning	2 values not saved	Write 'cas_settings.sav1'
iocdclv_S5R1									Warning	2 values not saved	Write 'cas_settings.sav1'
iocdclv_S5R2									Warning	2 values not saved	Write 'cas_settings.sav1'
iocdclv_S5R3									Warning	2 values not saved	Write 'cas_settings.sav1'
iocdclv_S6R1									Warning	2 values not saved	Write 'cas_settings.sav1'
iocdclv_S6R2									Warning	2 values not saved	Write 'cas_settings.sav1'
iocdclv_S6R3									Warning	2 values not saved	Write 'cas_settings.sav1'
iocgenFlasher	clonioc1.jlab.org	15 days, 00:19:01	1297141			10/25/2016 13:57:46			ok	ok	Write 'info_positions.sav1'
iocmoellerPS	clonioc1.jlab.org	15 days, 00:19:29	1297169			10/25/2016 13:57:18			ok	ok	Write 'info_positions.sav2'
iocmoellerTarget	clonioc1.jlab.org	14 days, 00:14:53	1210493			10/26/2016 14:01:54			ok	ok	Write 'info_positions.sav0'
iocmoellerCalib	clonioc1.jlab.org	7 days, 01:59:38	611978			11/02/2016 12:17:09			ok	ok	Write 'info_positions.sav1'
iocmoellerCalib	clonioc1.jlab.org	15 days, 00:18:59	1297139			10/25/2016 13:57:47			ok	ok	Write 'info_positions.sav1'
iocmoellerCalib	clonioc1.jlab.org	13 days, 01:00:28	1126828			10/27/2016 13:16:19			ok	ok	Write 'info_positions.sav1'
iocmoellerCalib	clonioc1.jlab.org	7 days, 00:16:44	605804			11/02/2016 14:00:02			ok	ok	Write 'info_positions.sav1'
iocmoellerCalib	clonioc1.jlab.org	15 days, 00:19:00	1297140			10/25/2016 13:57:47			ok	ok	Write 'info_positions.sav1'
iocgasSystem	clonioc1.jlab.org	7 days, 22:18:59	685139			11/01/2016 15:57:48			ok	ok	Write 'info_positions.sav1'
iocgasSystem86	svtsystem1.jlab.org	5 days, 02:40:39	441639			11/04/2016 11:36:08			ok	ok	Write 'info_positions.sav1'
iocagw	clonioc1.jlab.org	15 days, 00:19:00	1297140			10/25/2016 13:57:47			ok	ok	Write 'info_positions.sav1'
iocgenFlasher	clonioc2.jlab.org	43 days, 03:25:15	3727515			09/27/2016 10:51:32			Disconnect	Disconnected	Disconnected
iocmeCrates	clonioc2.jlab.org	43 days, 03:25:15	3727515			09/27/2016 10:51:32			Disconnect	Disconnected	Disconnected
iocmoellerPS	clonioc1.jlab.org	15 days, 00:19:00	1297140			10/25/2016 13:57:47			Disconnect	Disconnected	Disconnected
iocmoellerTarget	clonioc1.jlab.org	15 days, 00:19:00	1297140			10/25/2016 13:57:47			Disconnect	Disconnected	Disconnected
iocmoellerCalib	clonioc2.jlab.org	2 days, 23:05:37	255937			11/06/2016 14:11:09			Disconnect	Disconnected	Disconnected
iocmonitor	clonioc2.jlab.org	2 days, 23:05:37	255937			11/06/2016 14:11:09			Disconnect	Disconnected	Disconnected
VME IOCs											
IOC Name	Hostname	Up Time	Heartbeat	Expert	Reboot	Last Reboot	Status	Message	Autosave	Recently	Expert
iocclassc1	classc1	57 days, 02:12:12	4932732			09/13/2016 11:04:34	-	-	-	-	
iocclassc4	classc4	161 days, 17:16:03	13972563			05/31/2016 20:00:43	-	-	-	-	
iocclassc6	classc6	161 days, 17:14:06	13972447			05/31/2016 20:02:40	-	-	-	-	
iocclassc8	classc8	161 days, 17:14:06	13972447			05/31/2016 20:02:40	-	-	-	-	

FIG. 5: The primary IOC health screen showing uptime, heartbeats, and autosave status for each IOC, and buttons to restart them. Pink rows designate IOCs that are not currently running in this screenshot, but under normal operations all rows should have normal readings.

IV. High Voltage

The largest controls system in Hall B in terms of number of channels is high voltage (HV), with 23 CAEN crates including SY527, SY1527, and SY4527 mainframes. An overview screen of the status of all HV in Hall B is accessible from the HV button in the main CLAS12 menu as shown in Figure 6. Clicking on a detector in this overview screen will bring up the HV controls for that detector (also accessible under each detector's button in the main menu).

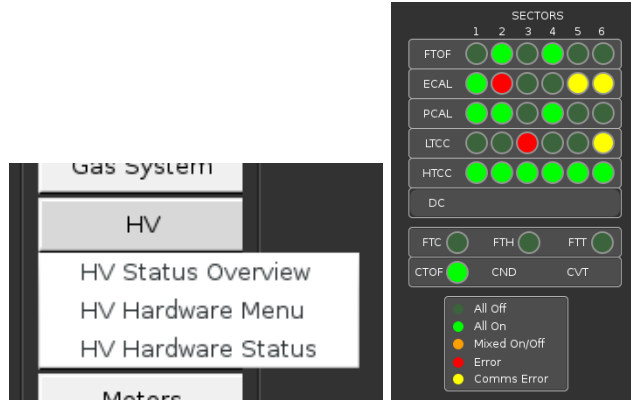


FIG. 6: Access to the HV overview screen from the main menu (left). Clicking on a detector's name in the overview screen (right) will open its HV controls screen.

V. Strip Charts

There are two applications available for plotting time histories of slow controls variables: StripTool and MyaViewer. Both are available from the *Strip Charts* button at the bottom of the main CLAS12 controls menu as shown in Figure 7.

The suggested tool for online operations in Hall B is StripTool, which has no access to archived data but is very robust and stable. MyaViewer is necessary for expert studies and can access the Mya archive used to store previous years of Hall B controls data. In either case, configuration files are loadable from their user interfaces to view a predetermined set of variables, or else you can choose any process variable to plot.

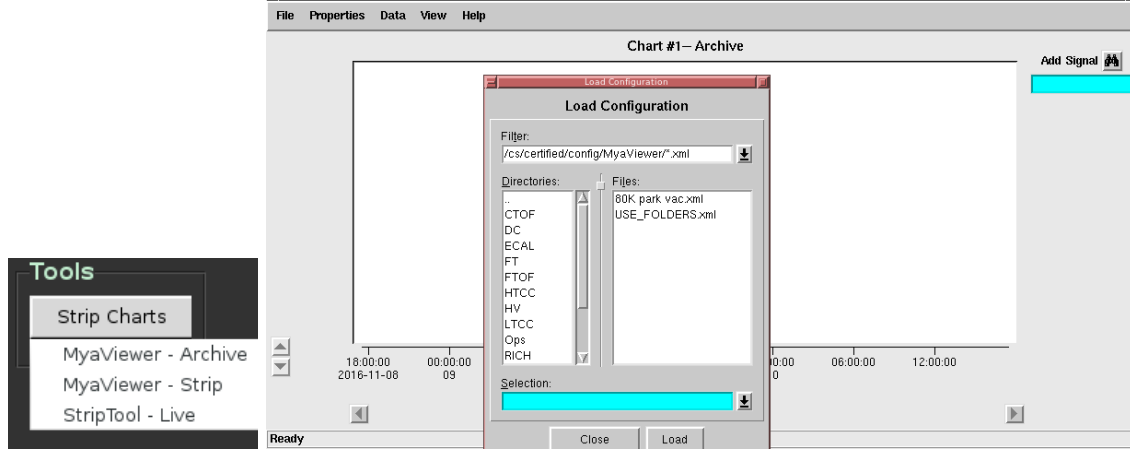


FIG. 7: Utilities for plotting time histories of slow controls variables are accessible from the *Tools* section of the CLAS12 main menu (left). An example of running MyaViewer and opening a preset configuration file via the *File* \rightarrow *Load Config* menu is shown on right.

VI. Logbook Entries and Screenshots

We use the JLab logbook system, and the primary Hall B logbook is called HBLOG and accessible in a web browser at

<https://logbooks.jlab.org/book/HBLOG>

In Hall B there are two primary methods for adding content to the logbook:

1. Use the web browser interface after logging in with your personal CUE credentials. That is the normal method used for filling out the shift checklist, updating a shift summary log entry, following up with comments on previous log entries, or adding more complex log entries (e.g. with multiple images).
2. Use our Hall B GUI that facilitates taking screenshots and quickly sending one to the HBLOG logbook as `user=clasrun`. This is accessed via the “logbook entry” item from the desktop menu, or via the following script in a terminal:

```
logbookEntry.sh
```

This is also the preferred method for taking screenshots and will always save them in `$HOME/screenshots` with timestamped filenames. See Figure 8 for details.

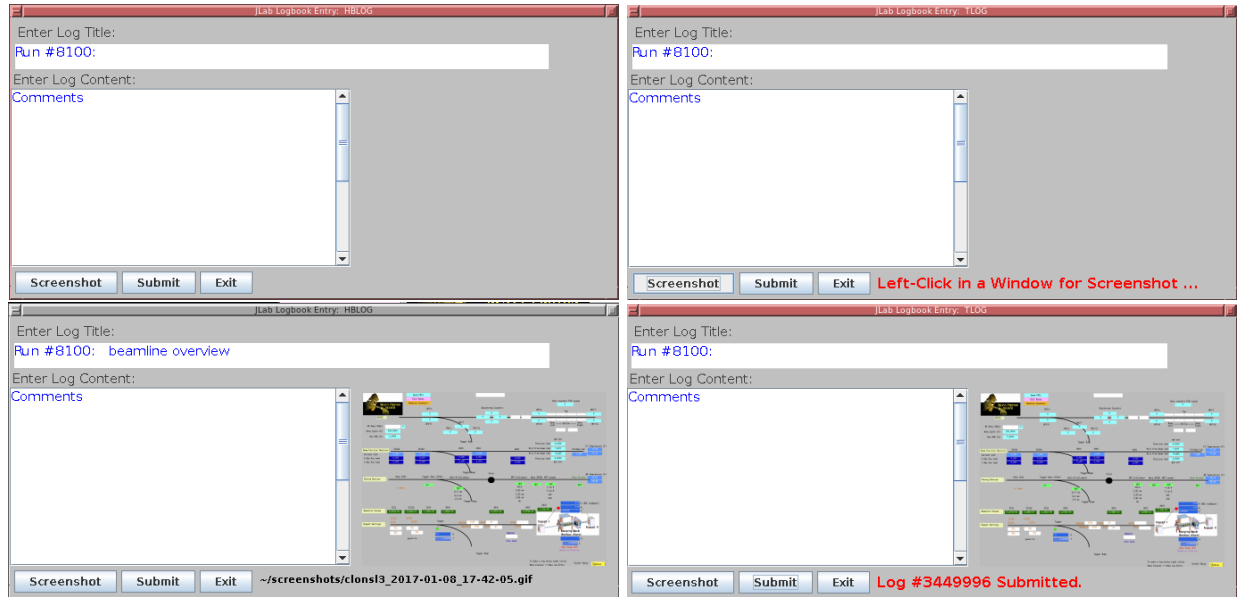


FIG. 8: Upon first opening the logbook/screenshot GUI (top left), only the log title has been automatically initialized (with the current run number). After clicking the “Screenshot” button (top right), it is waiting for you to left-click in the window you desire to capture (clicking the desktop instead of a window will capture the entire desktop). After taking a screenshot (bottom left), a snapshot of the image and its filename on disk are automatically displayed. Note that the “Screenshot” button can be used repeatedly to change the screenshot if you do not like the previous result, or just want to take more screenshots. The “Submit” button can be used to generate an entry in the HBLOG logbook, and after success the entry number will be displayed (bottom right).

VII. Accelerator Screens

The accelerator’s screens are accessed from the main CLAS12 menu via the *JMenu* button. This uses the `hbops` account on `h1b100`, a machine owned and maintained by the accelerator group. If a prompt requests a username, password, or terminal type, just press *Enter*. The location of the button on the CLAS12 menu and the *JMenu* screen that should appear are shown in Figure 9.

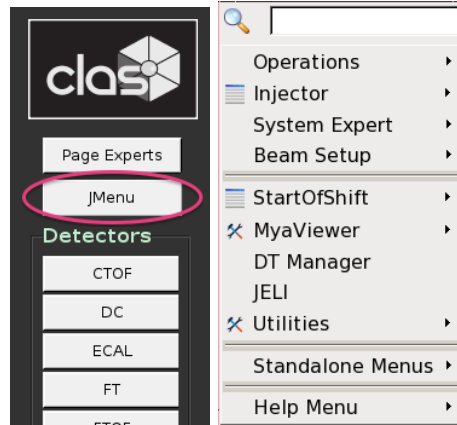


FIG. 9: The location of the button to access the accelerator screens from the CLAS12 controls menu (left) and the resulting accelerator JMenu main screen (right).

VIII. Paging System Experts

Paging on-call experts is available from the main CLAS12 controls menu via the *Page Experts* button at the very top of the screen. This will open a dropdown menu to choose the desired subsystem, and then open a new window in which to enter a message to be sent to the corresponding expert, as illustrated in Figure 10.

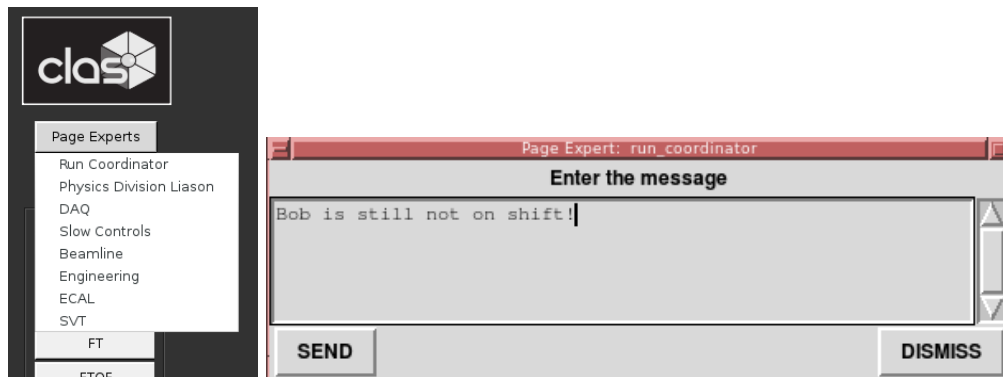


FIG. 10: The dropdown menu for choosing which expert to page (left) and the resulting dialog window in which to enter the message contents (right).

IX. Slow Controls Contacts

The individuals to be contacted for Hall B slow controls are shown in Table I. The first point of contact for shift operations is always the on-call controls expert, accessible from the paging system described in Section VIII of this document and the phone number in the first row of Table I. Additional contacts are listed in the table as a fallback.

On-Call		757-748-6922	
General	Nathan Baltzell	757-259-5902	baltzell@jlab.org
	Ken Livingston		kliv@jlab.org
	Wesley Moore	757-259-6033	wmoore@jlab.org
	Bryan McKinnon		mckinnon@jlab.org

TABLE I: Hall B slow controls contacts.

X. Remote Usage

There are separate server-grade machines for remote controls access, all with access to the same software and running the same operating system as the desktops in the counting house. For access outside the counting house, login to the server clonsl2. *In order to avoid heavily load on the machines used by counting house shift workers, it is important to not run on clonpc desktops remotely.* All controls computers are behind JLab's hallgw gateway and require 2-factor authentication for remote access.